

# ResMed Stellar 150 Download Instructions

Make sure the Stellar is on the home screen with the date and time.

**The humidifier has been removed in these pictures to simplify. It does not need to be removed in order to download.**



Remove the cap off the USB stick.

Insert the USB stick into either of the USB ports on the back of the machine

Now look at the screen for further instructions.



The USB stick will connect to the machine automatically in a few moments and on the screen will appear:

**“USB Data Management”**

Using the large black dial

Turn clockwise until the blue highlighted box is on the **“Write Set/Log Start”** box.

Push the black dial in to Select **“Write Set/Log Start”**




# ResMed Stellar 150 Download Instructions

The screen will change

Take the large black dial and turn clockwise until the blue highlighted box is on the “Detailed and Summary Data”

Push the black dial in to Select “Detailed and Summary Data”

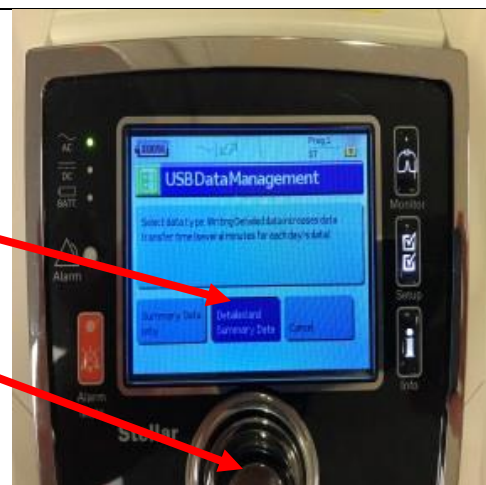

Turn the black dial to highlight “Select All”

Push the black dial in to accept. You will see green check marks beside all the dates.


Turn the black dial to highlight “Write”  
Push the black dial in to accept.



If the stick has previous data you may get a message “Overwrite Existing Data”

Take the black dial and turn to highlight “Yes” then  
Push the black dial in to accept.

Wait for the progress bar to complete  
(this may take several minutes)

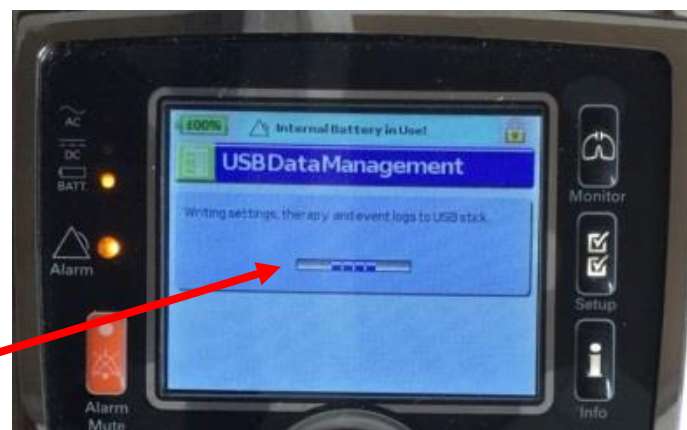


# ResMed Stellar 150 Download Instructions

When the download is complete a message will appear saying “written successfully to USB stick”.

In the bottom left corner “OK” is highlighted in blue.

Push the black dial in to accept.

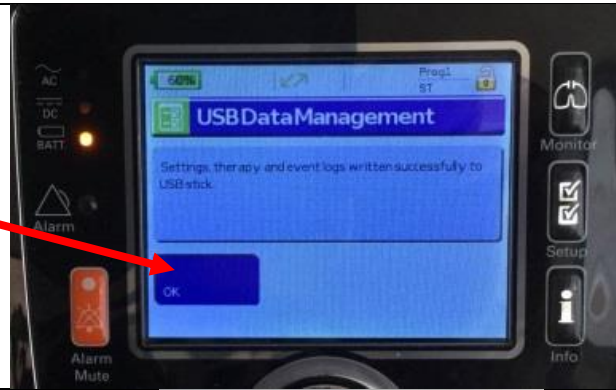

Remove the USB stick from the back of the unit.

Place the cap on the end of the USB stick.

Mail the USB stick to Dr. Leasa’s office.

